

## GRASP Visual Displays for Review

The GRASP experiment is designed to probe for changes in eye-hand coordination on orbit. Subjects (astronauts) perform a video game-like task using virtual reality that consists of shooting a projectile out of a hand-held tool toward a target. The target is presented as a set of parallel lines that can have different orientations in space and the tool and projectiles are presented as a line of spheres. Since the line of targets and projectiles are long and narrow, the tool must be placed at the same orientation as the target to achieve success.

GRASP studies the mental representation of the target's orientation by having subjects shoot the projectile to the **remembered** target orientation. The target is presented only briefly, then disappears. The subject then triggers the projectiles after a waiting period of about 5 seconds.

GRASP goes further to study how the brain represents objects in extra-corporal space. On most trials, the subject is required to move the head between the time that the target is present and the moment of the response. The subject must take into account that the object would no longer be at the same orientation with respect to the eyes when producing the response.

Trials vary in terms of what can seen and felt during the execution of the task, organized into three different experimental paradigms:

- In the main task (Visual-Manual) the subject can see the target but must trigger the projectiles without being able to see the orientation of the hand.
- In one variant (Visual-Visual) the subject can see both the target and the hand.
- In a second variant (Manual-Manual) the subject is not presented with an oriented visual target. Instead, he or she is guided by color cues to place the hand at the desired target orientation. The subject must memorize this orientation, lower the hand, wait for 5 seconds, and then reproduce the remembered orientation.

In practical terms, subjects perform repeated trials of one of these paradigms during each of three sub-sessions.

The subject also performs the experiment in two different postural conditions:

- On orbit, subjects are either strapped in an upright, seated position with respect to the ISS module, or for other trials, are maintained in a state of quasi-freefloating with no

contact with the floor, walls or ceiling of the ISS module.

- On the ground subjects perform the experiment seated upright in a chair or lying supine with respect to gravity.

The details of the experiment design can be found elsewhere (e.g. the ESR). In this document we present the different visual displays that subjects and operators will see when performing the GRASP experiment. The document is organized as follows:


- This Introduction
- Graphical User Interface that controls the experiment on the Perspectives laptop computer.
- Instructions presented to the operator during hardware preparation.
- Screens for the Visual-Visual Task
  - Instructions (presented in the GUI)
  - VR Scenes (presented in the Virtual Reality Helmet)
- Screens for the Manual-Manual Task
  - Instructions (presented in the GUI)
  - VR Scenes (presented in the Virtual Reality Helmet)
- Screens for the Visual-Manual Task
  - Instructions (presented in the GUI)
  - VR Scenes (presented in the Virtual Reality Helmet)
- Additional messages that the subject may see inside the VR display.

# GRASP Software User Interface


This document provides an overview of the steps followed by a subject or operator in order to run the GRASP installation and scientific protocols on the Perspectives laptop. It covers only the software aspects, assuming that the GRASP hardware (including shared GRIP hardware) is already installed per the PODF specific to those activities.

## ***Starting GRASP on the Perspectives Laptop***

A crew member initiates Grasp activities by double clicking on the GRASP icon situated on the desktop. This will start up first a utility program to verify the installation of the virtual reality engine and VR Headset. The initial screen will be as follows while the VR system initializes:

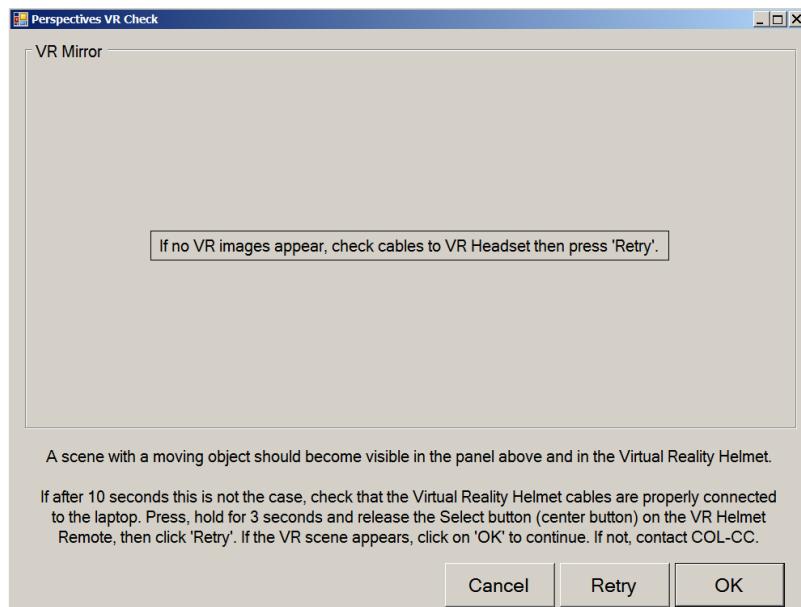


after 10-15 seconds a VR screen should appear as follows:

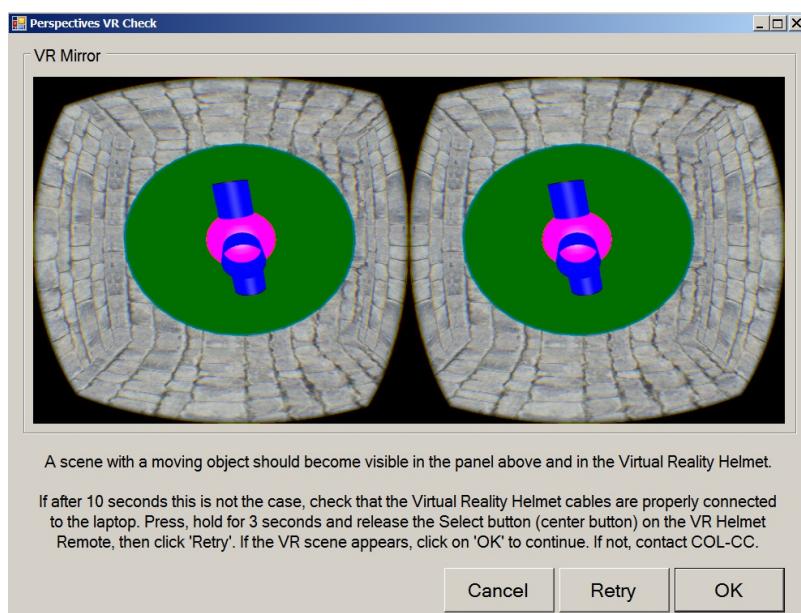


If the VR scene does not appear after 30 seconds, the operator should press and hold for 3 seconds the Select button on the VR Remote and then release, as indicated by the instructions on the screen. This should cause the VR scene to appear in the window and in the VR Headset. If not, the operator should contact Col-CC for further instructions.

After the VR engine has been started successfully, a window showing the status of the 3D tracker:

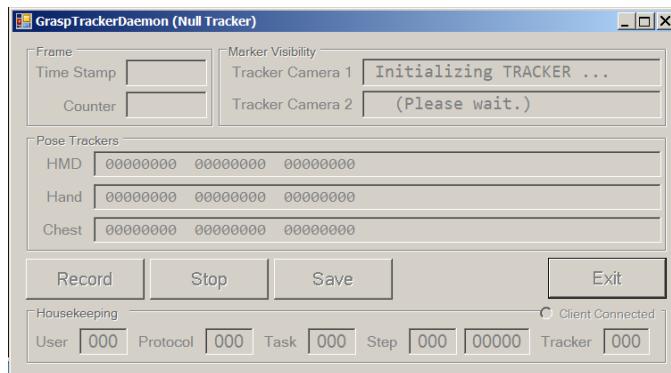


should appear for ~10 sec and then disappear (collapse to task bar).

Finally, the GRASP main menu should appear:

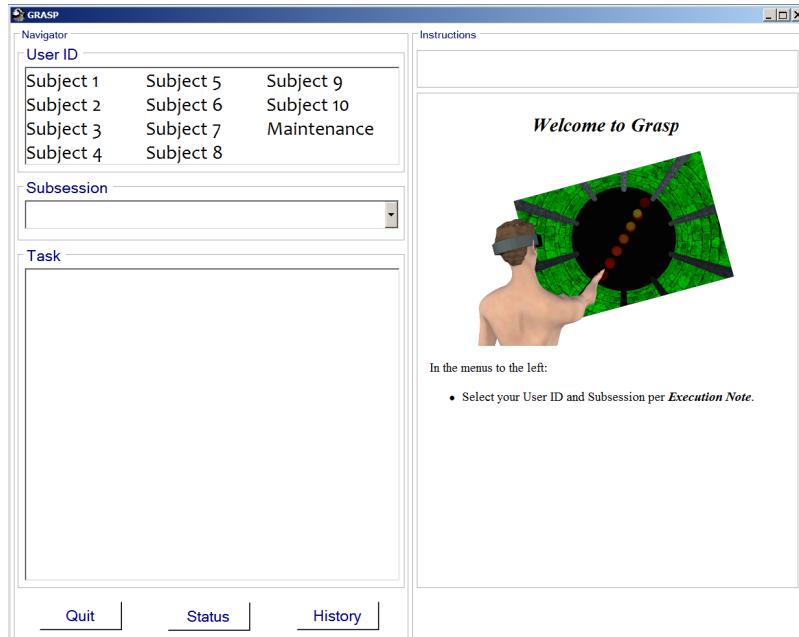


The subject or operator selects his or her designated subject ID in the upper-left panel of GRASP@ISS. Then the desired sub-session ID sub-session is selected via a pull-down menu:

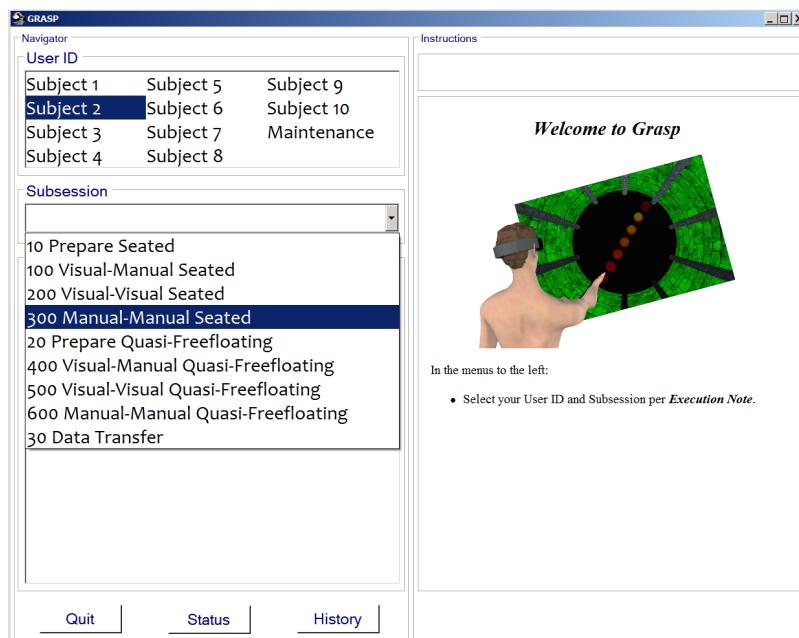


The subject ID and the sub-session to be performed at any given time is communicated to the crewmember via an Execution Note.

Selecting the sub-session protocol brings up a list of tasks to be executed in a defined order. The first task is automatically selected. An example is shown here:

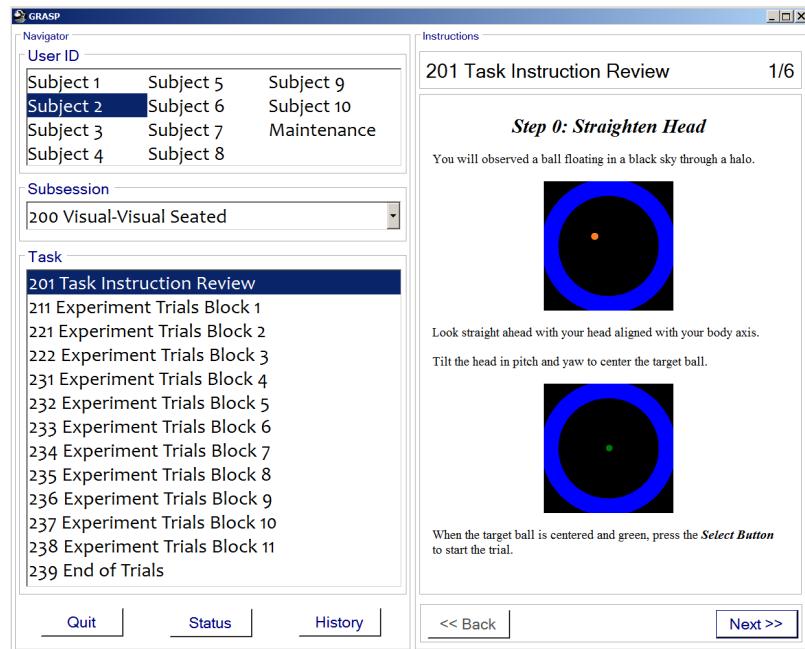

Note that at any time the subject can press the 'Status' button to bring up a screen showing the visibility of the markers on various marker structures and a test window showing the status of the VR hardware.

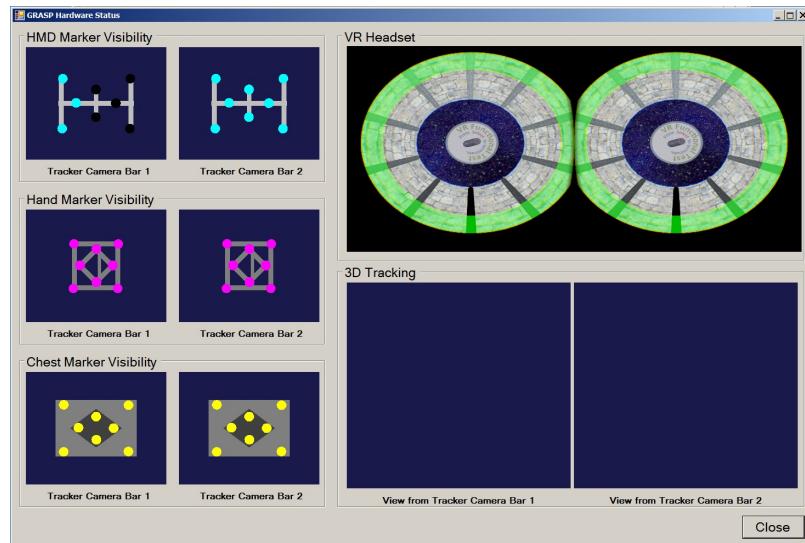

The subject can also press the 'History' button to display a list of the tasks that have been carried out on that day, with start and end time and an indication of the completion code for each task.

Start	Usr	Ssn	Tsk	Stp	Completion	Status	End
2017.10.16 00:37:11	099	950	952	001	Completed normally	(code 0)	2017.10.16 00:37:11
2017.10.16 00:39:06	099	950	952	001	Completed normally	(code 0)	2017.10.16 00:39:06
2017.10.16 00:40:21	099	950	952	001	Completed normally	(code 0)	2017.10.16 00:40:22
2017.10.16 00:41:03	099	950	952	001	Completed normally	(code 0)	2017.10.16 00:41:03
2017.10.16 00:41:43	099	950	952	001	Completed normally	(code 0)	2017.10.16 00:41:43
2017.10.16 00:46:03	099	950	951	001	Completed normally	(code 0)	2017.10.16 00:46:03
2017.10.16 00:47:11	099	950	952	001	Completed normally	(code 0)	2017.10.16 00:47:11
2017.10.16 01:31:42	099	950	954	001	Completed normally	(code 0)	2017.10.16 01:31:42
2017.10.16 01:32:07	099	950	955	001	Completed normally	(code 0)	2017.10.16 01:32:07
2017.10.16 01:34:16	099	950	954	001	Completed normally	(code 0)	2017.10.16 01:34:16
2017.10.16 01:34:56	099	950	955	001	Completed normally	(code 0)	2017.10.16 01:34:57
2017.10.16 01:35:44	099	950	952	001	Completed normally	(code 0)	2017.10.16 01:35:44

## Executing Tasks

The panel on the right then guides the crew member through the steps required by each task. The panel may take one of four different forms.


1. First, the crew member is presented with one or more pages of instructions. Once the instructions have been read the subject proceeds to the subsequent page by pressing the 'Next' button. It is also possible to return to the previous step (except when at the first step) by pressing the 'Back' button.

Instructions


### 202 Task Instruction Review 1/5

**Step 1: Straighten Head**

- If your head is laterally tilted with respect to the trunk, the halo will be red.

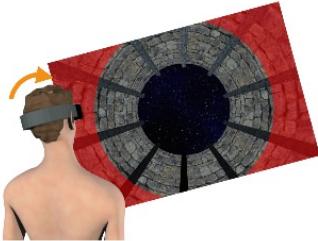

- Align the head to your body axis, this will make the halo turn bright green.

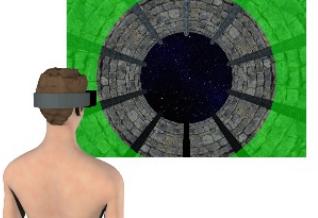

<< Back Next >>

2. The crew member may then be prompted to launch a particular action. In the example shown here, the system is ready to perform the tracker alignment. The crew member presses the 'Execute' button to initiate the action. He or she may also press 'Back' to return to the previous action or instruction (except when it is the first step in a task). In exceptional cases, the subject may be told by COL-CC to press 'Skip' to move on to the next step without performing the action.


Instructions

11 Preparation 1/4

## Tracker Alignment

Prepare to perform the alignment procedure for the 3D tracker:

- Check that the **GRIP Chair** is installed with the seat back deployed and the **Chest Marker Support**, **Hand Marker Support** and **Head Marker Support** are attached with velcro to the **GRIP Chair** as shown.

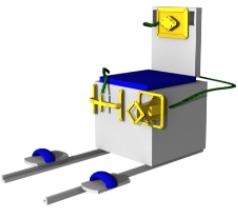

**Note:** Alignment of the **Chest Marker Support** on the backrest is critical. Please ensure that the top edge of the **Chest Marker Support** is aligned with the top edge of the backrest with the cable exiting to the right side of the chair (subject's right side when seated facing the **Tracking Cameras**).

- Check that the field of view of the **Tracking Cameras** is free of obstructions.

Press '**Start**' to perform alignment.

<< Back Skip >> Start

3. Pressing 'Execute' typically causes a new program to execute, which will activate specific popup windows outside the GRASP@ISS main window. When execution of the task terminates, focus will return to the GRASP@ISS main window and the termination status of the task is displayed. If the task was successful, the crew member will typically see the following message allowing them to confirm and move on to the following step (by pressing 'Next') or to return and repeat the task (by pressing 'Back').



4. If the task did not complete successfully, an alert will be displayed that shows the error return code (in the 'Error Code' box) and that typically includes a graphical or text explanation of the anomaly. Shown here is the generic message that may be used for any anomalous condition, but more specific messages may also be displayed (see below).

Instructions

203 Practice Trials 1/1

**Action Completed with Errors**



Press:

- '**Retry**' to repeat this action.
- '**Ignore**' to ignore and move on.
- '**Restart**' to restart this task or select another.

**Recommendation:** Attempt '**Retry**' no more than once or twice before consulting with POCC. '**Ignore**' and '**Restart**' should not be used without POCC authorization.

Error Code  **Restart** **Ignore** **Retry**


In these circumstances, the crew member is invited to repeat the task by pressing 'Retry'. If a second attempt also fails, the crew member has the option to move on to the next step anyway, by pressing 'Ignore' or to repeat the entire task from the begining by pressing 'Restart'.

The full set of possible instruction, status and alert screens are included in the following sections.

## Tracker Alignment

Prepare to perform the alignment procedure for the 3D tracker:

- Check that the ***GRIP Chair*** is installed with the seat back deployed and the ***Chest Marker Support***, ***Hand Marker Support*** and ***Head Marker Support*** are attached with velcro to the ***GRIP Chair*** as shown.

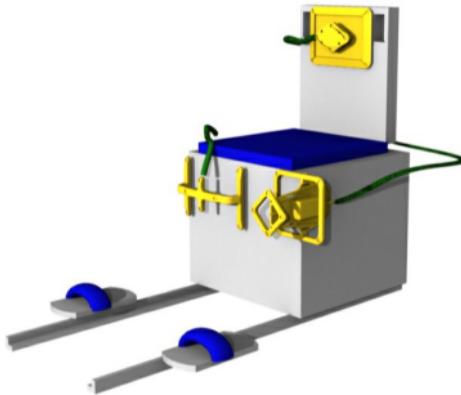


**Note:** Alignment of the ***Chest Marker Support*** on the backrest is critical. Please ensure that the top edge of the ***Chest Marker Support*** is aligned with the top edge of the backrest with the cable exiting to the right side of the chair (subject's right side when seated facing the ***Tracking Cameras***).

- Check that the field of view of the ***Tracking Cameras*** is free of obstructions.

Press '**Start**' to perform alignment.


# Tracker Alignment Running

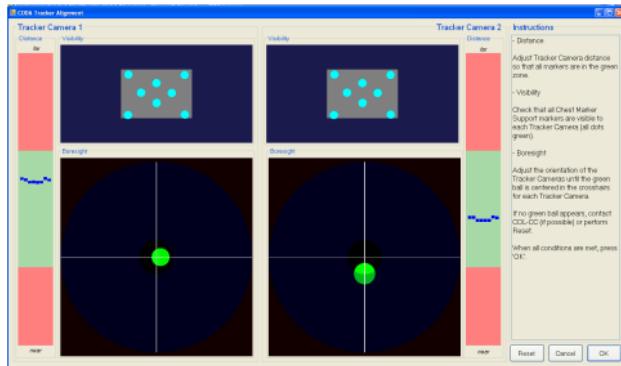

**Follow the instructions on the screen.**

## HMD Assembly

- Attach the ***Head Marker Support*** to the front of the ***Virtual Reality Headset*** as shown:

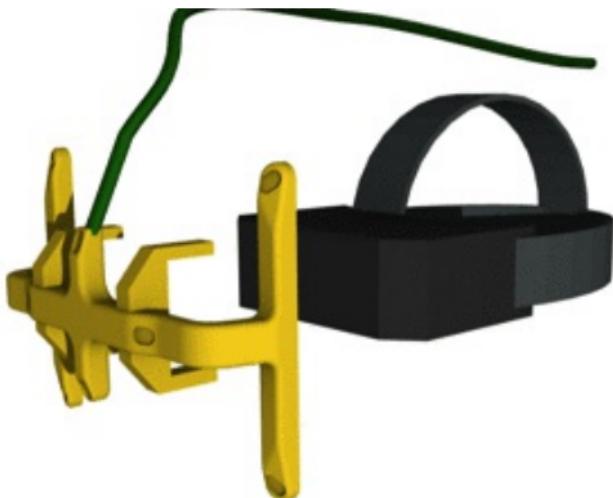

- Secure with velcro straps on top and bottom.
- Attach the small junction box on the ***Head Marker Support*** cable to the top strap of the ***Virtual Reality Headset***.


**Note:** Alignment of the Marker Structure to the ***Virtual Reality Headset*** is critical. Please ensure that the structure is parallel to the front face of the ***Virtual Reality Headset***.

Press '**Next**' to continue.

## Objects Pre-Positioning

Place the objects as shown to prepare for experiment execution:

- Replace the Cushion on the ***GRIP Chair***.
- Temp stow the ***Virtual Reality Headset***, the ***Hand Marker Support*** and the ***Chest Marker Support*** to the right side of the ***GRIP Chair*** (subject's right side when seated facing the ***Tracking Cameras***).

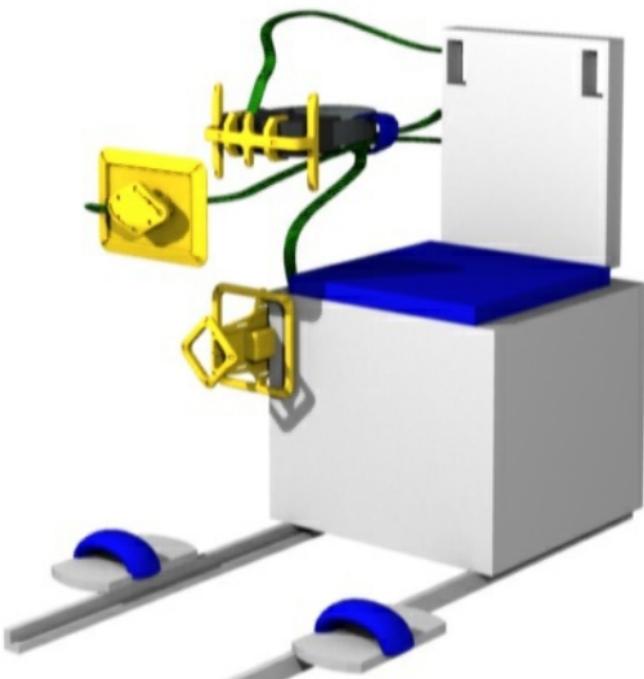


- Ensure that the ***Marker Array Cables*** and ***VR Headset Cable*** pass behind the ***GRIP Chair***.

Press '**Next**' to continue.

# Tracker Alignment

Prepare to perform the alignment procedure for the 3D tracker:

- Check that the quasi-freefloating ***Floating Backrest*** is installed with the ***Chest Marker Support, Hand Marker Support*** and ***Head Marker Support*** attached with velcro as shown to the restraint plate.

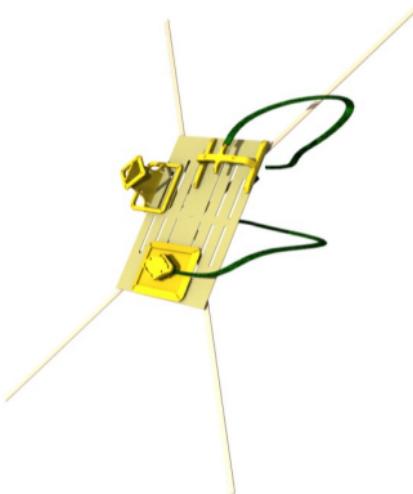


**Note:** Alignment of the ***Chest Marker Support*** on the restraint plate is critical. Check that the top edge of the ***Chest Marker Support*** is aligned with the edge of the plate nearest the deck of the ISS module and that the cable exits the support to your right when facing the plate.

- Check that the field of view of the ***Tracking Cameras*** is free of obstructions.

Press '**Start**' to perform alignment.


# Tracker Alignment Running

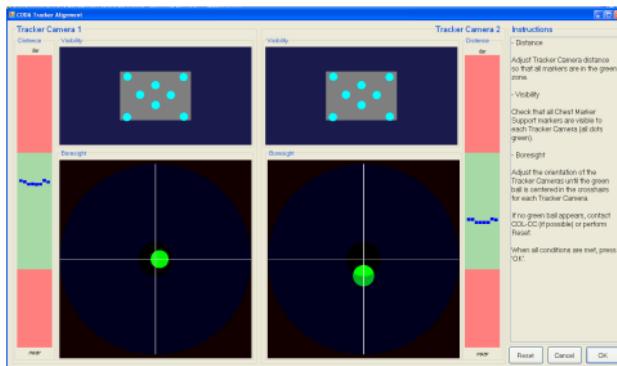

**Follow the instructions on the screen.**

## HMD Assembly

- Attach the ***Head Marker Support*** to the front of the ***Virtual Reality Headset*** as shown:

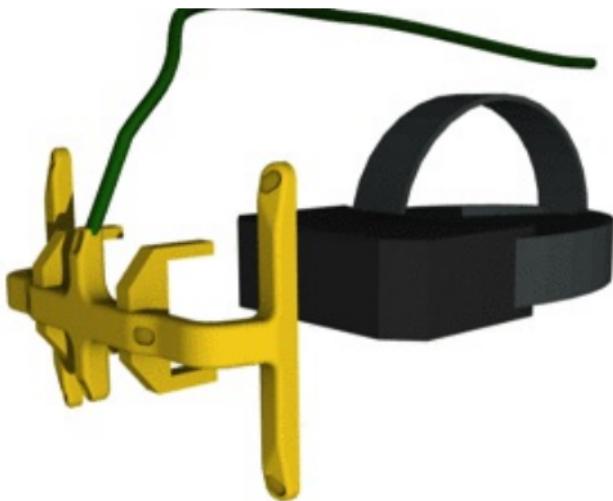

- Secure with velcro straps on top and bottom.
- Attach the small junction box on the ***Head Marker Support*** cable to the top strap of the ***Virtual Reality Headset***.


**Note:** Alignment of the Marker Structure to the ***Virtual Reality Headset*** is critical. Please ensure that the structure is parallel to the front face of the ***Virtual Reality Headset***.

Press '**Next**' to continue.

## Objects Pre-Positioning

Place the objects as shown to prepare for experiment execution:

- Replace the Cushion on the ***GRIP Chair***.
- Temp stow the ***Virtual Reality Headset***, the ***Hand Marker Support*** and the ***Chest Marker Support*** to the right side of the ***GRIP Chair*** (subject's right side when seated facing the ***Tracking Cameras***).

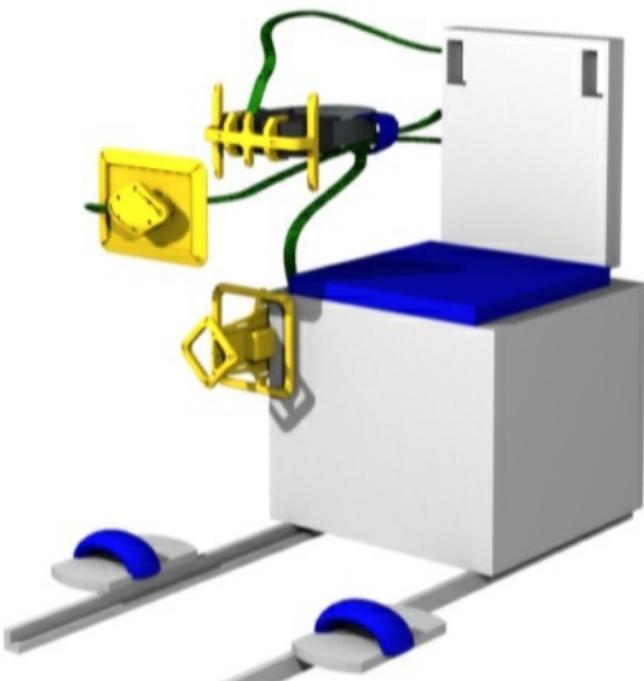

- Ensure that the ***Marker Array Cables*** and ***VR Headset Cable*** pass behind the ***GRIP Chair***.

Press 'Next' to continue.


## Visual-Visual

You have selected the Visual-Visual task.


Press '**Next**' to receive further instructions.

## Step 0: Straighten Head on Shoulders

You will observe a ball floating in a black sky through a halo.

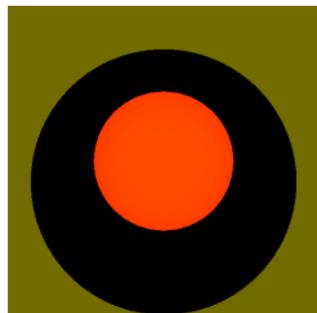

Look straight ahead with your head aligned with your body axis.  
Tilt your head in pitch and yaw until the ball turns green.  
Tilt your head in roll until the halo turns bright green.

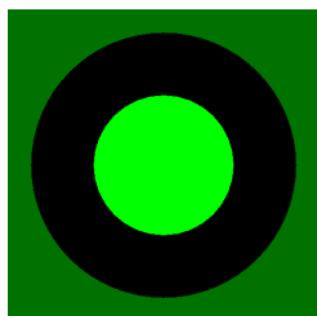


When the ball and halo are green, press the **Select Button** to start.




## Step 1: Straighten Head


If the halo is red you must reorient your head:

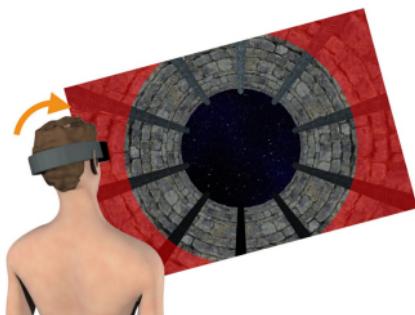

Look straight ahead then slowly tilt your head on your shoulders. This will make the halo turn progressively from red to green:

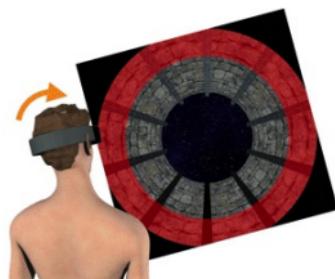


When the halo pops to bright green, hold this position.

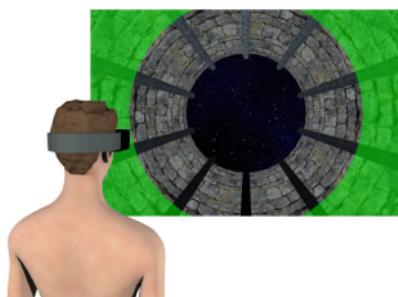

## ***Step 2: Acquire the Target Orientation***

Once your head is aligned with your body (i.e. when the halo is bright green), a set of parallel bars will appear at the end of the tunnel, representing the target orientation.

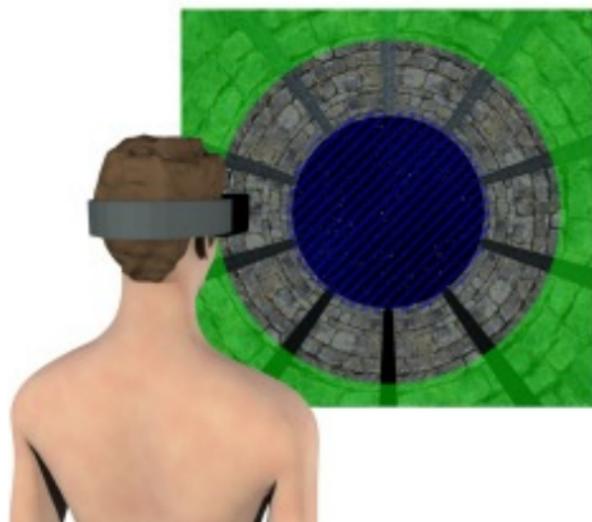


Memorise the target orientation.


**Note:** Your hand must remain in the neutral position on your lap throughout the target presentation. Please ***do not*** align your hand to the target orientation during target presentation or during the memory delay period.

### *Step 3: Tilt the Head*

If the halo turns red, you must tilt your head to a new orientation. The arrow tells you which way to rotate the head.

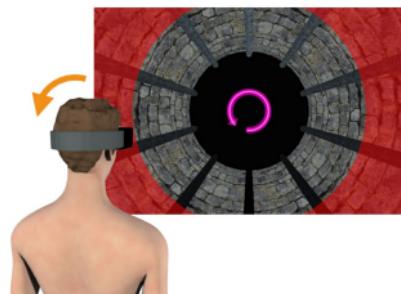


Roll your head slowly in the indicated direction until the halo turns bright green and hold this position.

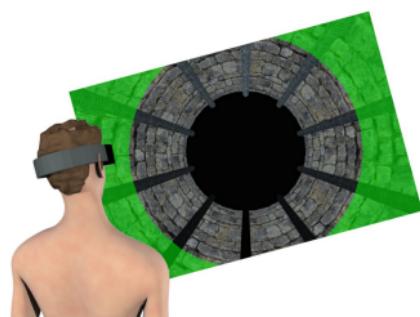

## Step 4: Align to Remembered Target

With the left hand use the **VR Headset Remote** to align the virtual tool to the memorised target.

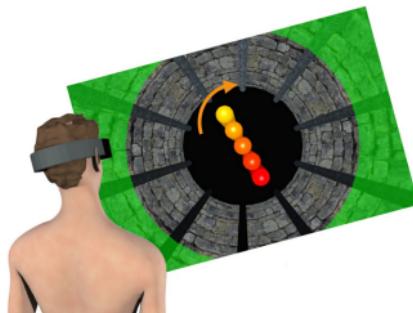

**Note:** Press and hold RIGHT and LEFT on the **VR Headset Remote** to rotate the virtual tool clockwise and counter-clockwise, respectively.




Press the **Select Button** on the **VR Headset Remote** to validate your response once you think that the virtual tool is aligned with the memorised target.


## ***Step 5: Check Response***

In some of the trials a line of balls will be projected from your hand towards the targets, allowing you to verify whether or not you correctly aligned your hand. On other trials you will see only a single spherical projectile without any indication as to whether your hand was correctly aligned to the memorised target orientation.


A green sphere appears briefly if the task sequence is completed correctly, regardless of the accuracy of the hand's alignment to the target.

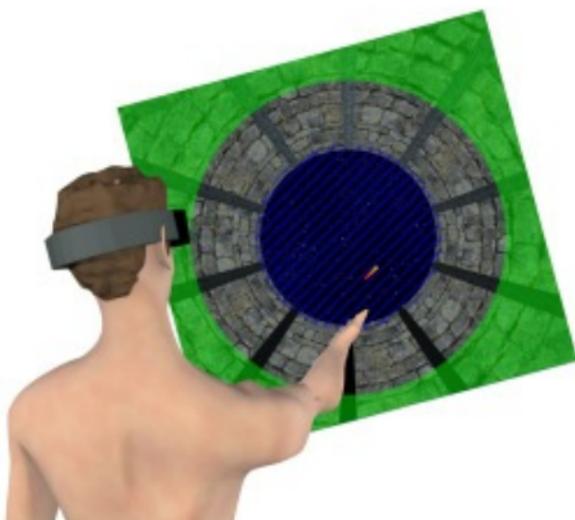

The Visual-Visual task starts with a welcome screen that prompts the subject to press a button in order to start a block of trials:

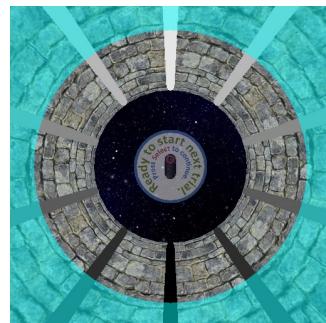


The subject is prompted to look straight ahead and bring the head to an upright position with respect to the shoulders.

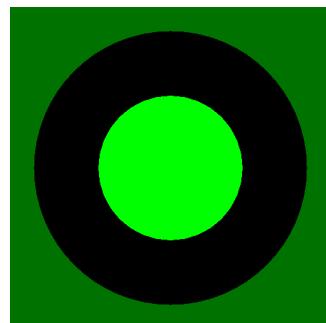


The subject should straighten the head first, then press the Select button on the VR Remote to continue.

A target orientation is then presented to the subject in the form of grill of parallel lines situated at the end of the tunnel:

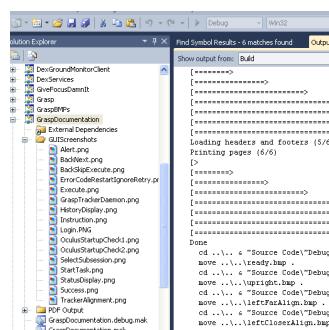

The target disappears and the subject is prompted to tilt the head to a new orientation through the same system of halo colors and direction arrows as described previously:

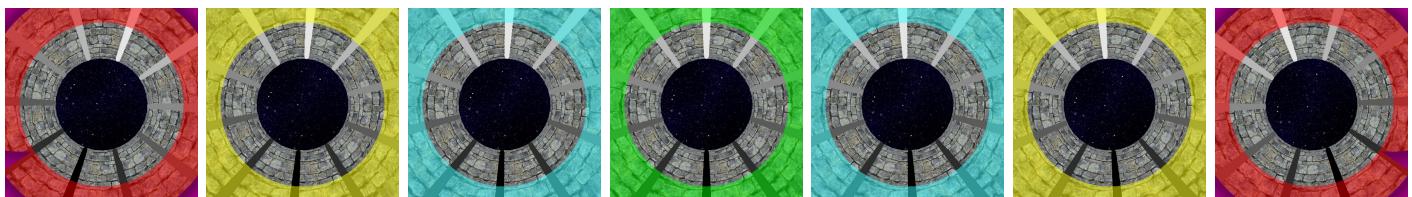


Tilted to the Left ----- Correct Orientation ----- Tilted to the Right

The subject then reproduces the remembered target orientation by turning the virtual tool using the Remote Control. The subject presses the right and left arrows of the Navigation Ring on the Remote Control to turn the virtual tool clockwise or counter-clockwise, respectively.




The subject presses the Select Button on the VR Remote to validate the trial and move on to the next. If the block of trials is completed, a message is displayed to let the subject know:

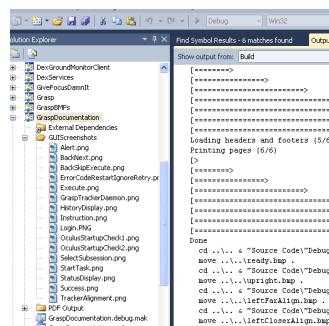

The subject may press the Select button on the VR Remote to exit the VR scene and pause the experiment. Or he or she can press and hold the Return button on the VR Remote to continue to the next block of trials.

If all blocks in the subsession are completed, the subject is prompted to doff the Virtual Reality Headset.

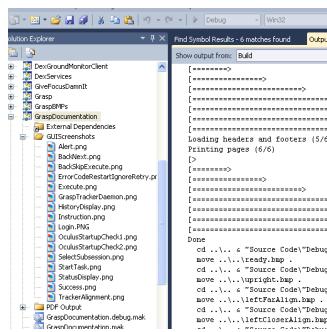

Note that other circular prompt with text instructions may appear at different times in response to different conditions or actions by the subject. For instance, a message will be generated if the subject does not maintain the head in the specified tilted position. All possible messages are shown in a separate section of this document.


## **Manual-Manual**

You have selected the Manual-Manual task.


Press '**Next**' to receive further instructions.

## Step 0: Straighten Head on Shoulders

You will observe a ball floating in a black sky through a halo.

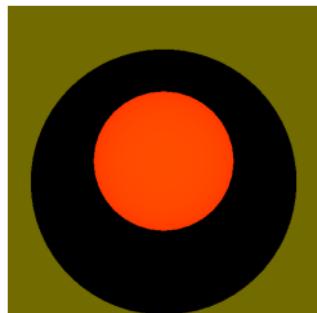

Look straight ahead with your head aligned with your body axis.  
Tilt your head in pitch and yaw until the ball turns green.  
Tilt your head in roll until the halo turns bright green.

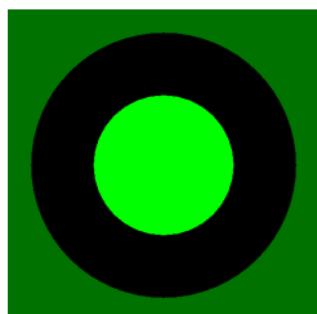


When the ball and halo are green, press the **Select Button** to start.




## Step 1: Straighten Head


If the halo is red you must reorient your head:

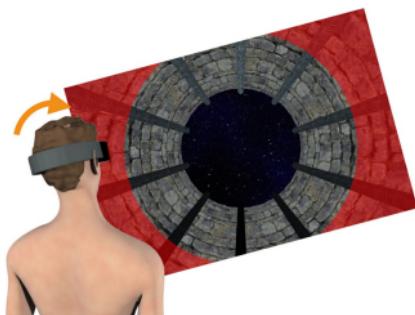

Look straight ahead then slowly tilt your head on your shoulders. This will make the halo turn progressively from red to green:

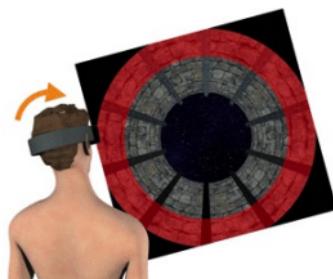


When the halo pops to bright green, hold this position.

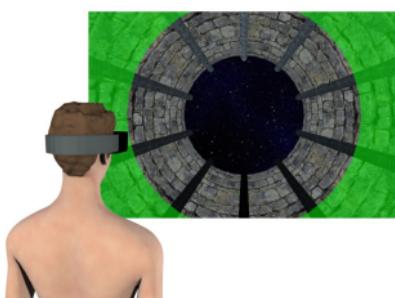


## Step 2: Acquire the Target Orientation

Once your head is aligned with your body (i.e. when the halo is bright green), a yellow spherical target will appear. Raise your arm and point your hand such that the virtual laser-pointer falls on the target.

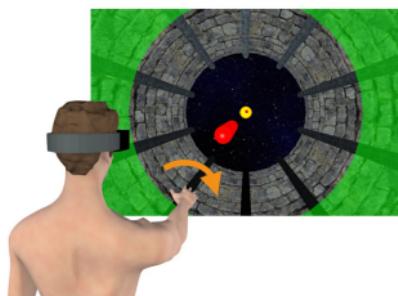

Rotate your hand in a rolling motion around the axis of your arm until the virtual pointer and tool turn bright green. Memorise the hand orientation.

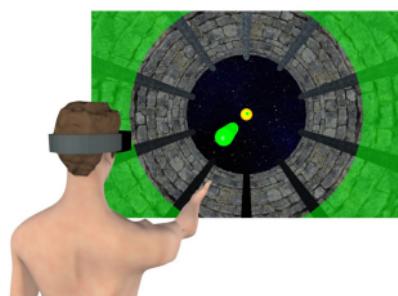


When the target disappears, lower your arm to a neutral position on your lap.


**Note:** Please **do not** try to maintain the memorised roll orientation of the hand. After target acquisition you should return your hand to the same neutral position on every trial.

### *Step 3: Tilt the Head*

If the halo turns red, you must tilt your head to a new orientation. The arrow tells you which way to rotate the head.

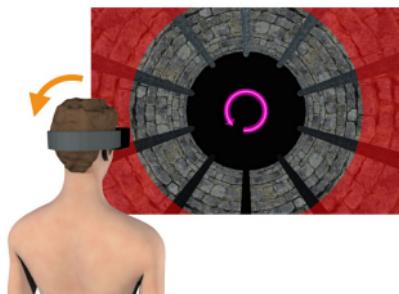


Roll your head slowly in the indicated direction until the halo turns bright green and hold this position.

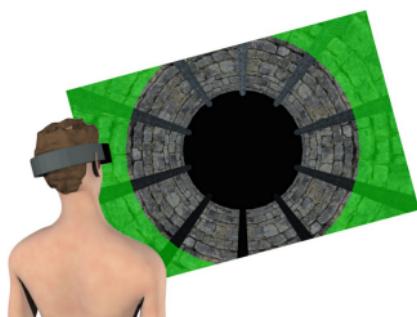

## Step 4: Align to Target

When the orange disk reappears, raise your arm and aim with the hand until the virtual laser pointer falls on the disk. Rotate your hand in a rolling motion around the axis of your arm to the remembered target orientation.

**Note:** In this phase the virtual tool gives no indication about the orientation of the hand around the roll axis. You must rotate the hand to the **remembered** orientation while pointing to the target in azimuth and elevation.

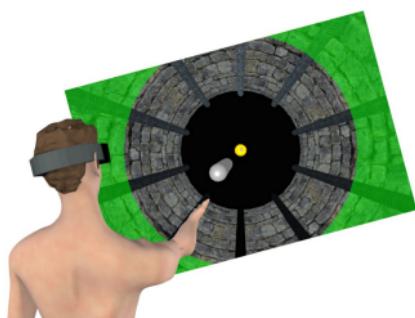


With your left hand, press the **Select Button** on the **VR Headset Remote** to validate your response once you think your hand is aligned with the memorised target.




## *Step 5: Check Response*

In some of the trials a line of balls will be projected from your hand towards the targets, allowing you to verify whether or not you correctly aligned your hand. On other trials you will see only a single spherical projectile without any indication as to whether your hand was correctly aligned to the memorised target orientation.


A green sphere appears briefly if the task sequence is completed correctly, regardless of the accuracy of the hand's alignment to the target.

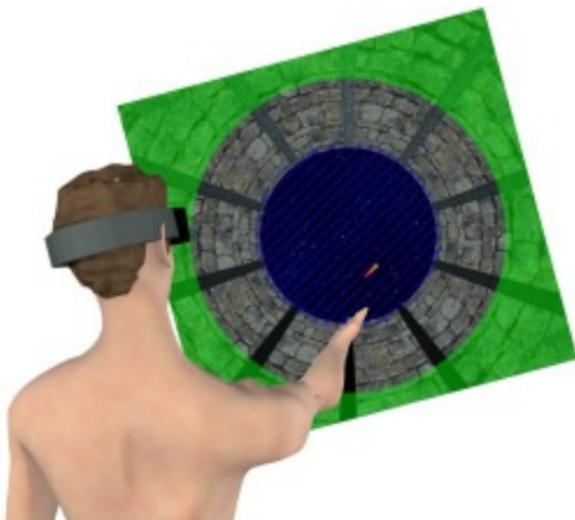

The Manual-Manual task starts with a welcome screen that prompts the subject to press a button in order to start a block of trials:

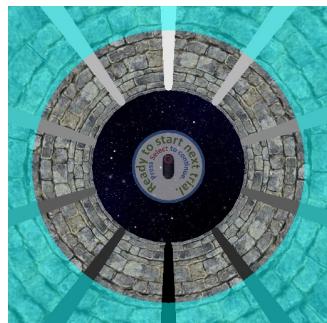


The subject is prompted to look straight ahead and bring the head to an upright position with respect to the shoulders.

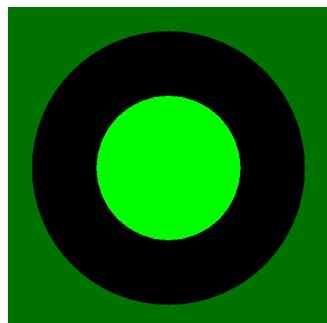


The subject should straighten the head first, the press the Select button on the VR Remote to continue.

The subject is then guided to place the hand at the target orientation. The subject sees a) an orange spherical target at the end of the tunnel, b) a virtual representation of a tool in the visual scene (a cylinder) that moves with the hand in 3D visual space and c) a 'laser pointer' that projects from the tool along the axis of the fingers. The virtual tool changes color according to the roll orientation of the hand (pronation/supination). The subject must raise the arm in front of the body, aim the tool by moving the arm in azimuth and elevation until the laser pointer falls on the spherical target and then rotate the hand until the tool turns green. If the orientation of the hand is far from the target orientation, an arrow appears around the virtual tool showing which way to rotate the hand.

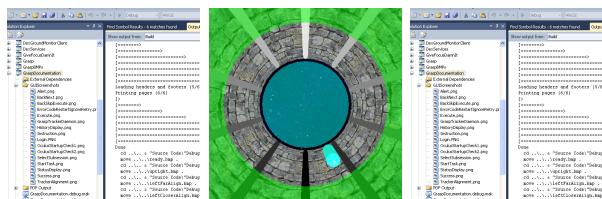


Once the target orientation has been achieved and held for 1 second, the tool and laser pointer disappear. The subject must lower the arm and is then prompted to tilt the head to a new orientation through the same system of halo colors and direction arrows as described previously:

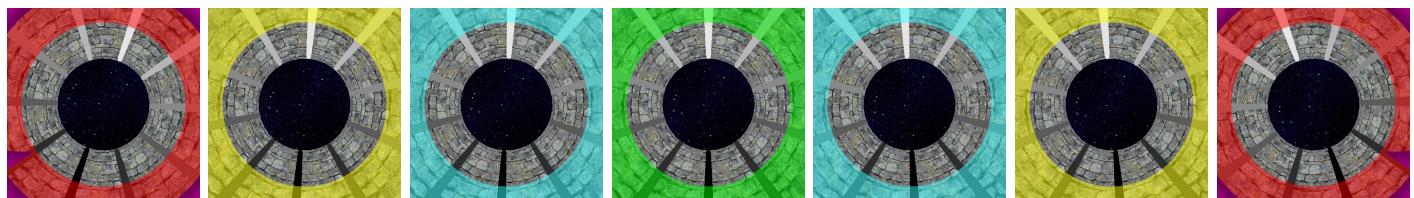


Tilted to the Left ----- Correct Orientation ----- Tilted to the Right

Finally subject then reproduces the remembered target orientation with the hand. The virtual tool, laser pointer and spherical target reappear, but this time the tool and laser pointer are constantly blue. The subject must raise the hand and aim with the laser pointer at the spherical target but then rotate the hand to the remembered orientation without any further visual cues:

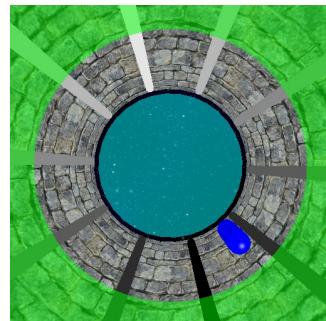


The subject presses the Select Button on the VR Remote to validate the trial and move on to the next. If the block of trials is completed, a message is displayed to let the subject know:

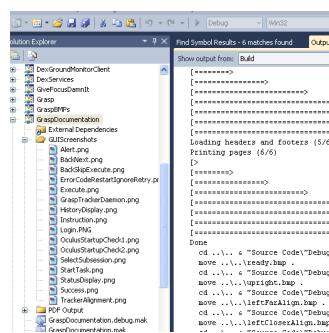

The subject may press the Select button on the VR Remote to exit the VR scene and pause the experiment. Or he or she can press and hold the Return button on the VR Remote to

continue to the next block of trials.

If all blocks in the subsession are completed, the subject is prompted to doff the Virtual Reality Headset.

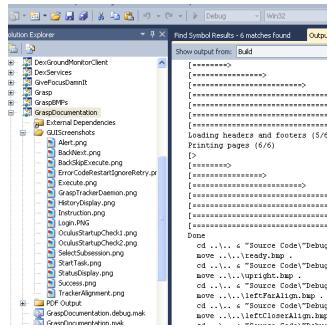

Note that other circular prompt with text instructions may appear at different times in response to different conditions or actions by the subject. For instance, a message will be generated if the subject does not maintain the head in the specified tilted position. All possible messages are shown in a separate section of this document.


# Visual-Manual

You have selected the Visual-Manual task.


Press '**Next**' to receive further instructions.

## Step 0: Straighten Head on Shoulders

You will observe a ball floating in a black sky through a halo.

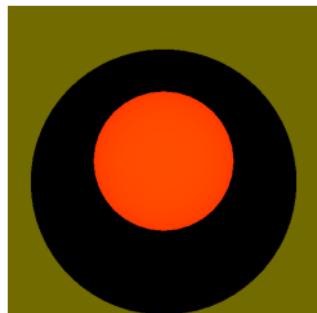

Look straight ahead with your head aligned with your body axis.  
Tilt your head in pitch and yaw until the ball turns green.  
Tilt your head in roll until the halo turns bright green.

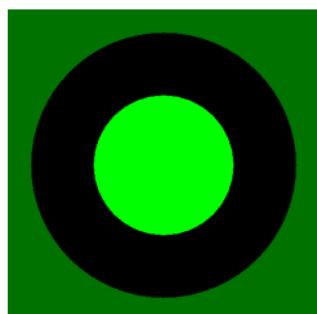


When the ball and halo are green, press the **Select Button** to start.




## Step 1: Straighten Head


If the halo is red you must reorient your head:

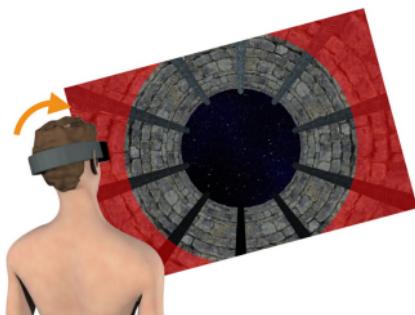

Look straight ahead then slowly tilt your head on your shoulders. This will make the halo turn progressively from red to green:

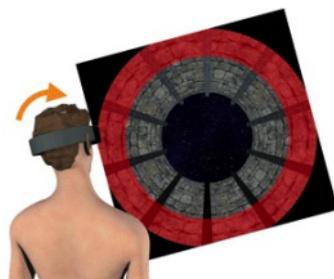


When the halo pops to bright green, hold this position.

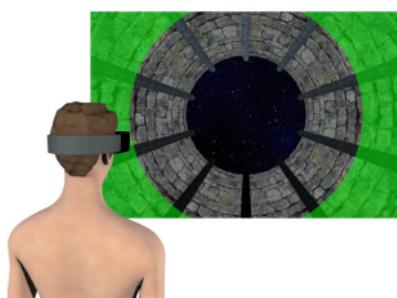

## *Step 2: Acquire the Target Orientation*

Once your head is aligned with your body (i.e. when the halo is bright green), a set of parallel bars will appear at the end of the tunnel, representing the target orientation.

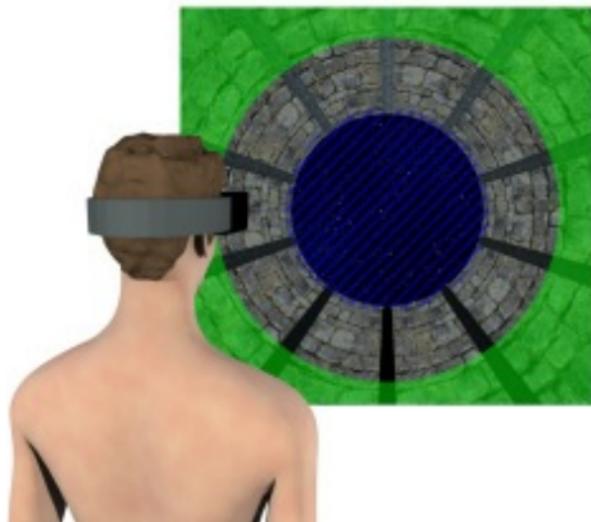


Memorise the target orientation.


**Note:** Your hand must remain in the neutral position on your lap throughout the target presentation. Please ***do not*** align your hand to the target orientation during target presentation or during the memory delay period.

## *Step 3: Tilt the Head*

If the halo turns red, you must tilt your head to a new orientation. The arrow tells you which way to rotate the head.

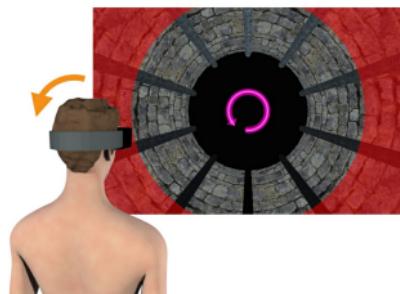


Roll your head slowly in the indicated direction until the halo turns bright green and hold this position.

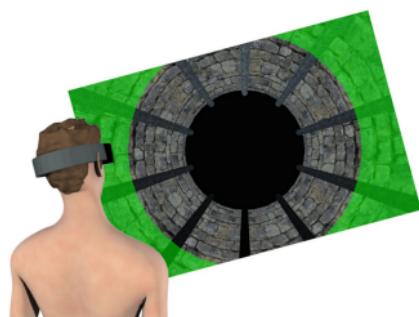

## Step 4: Align to Target

When the orange disk reappears, raise your arm and aim with the hand until the virtual laser pointer falls on the disk. Rotate your hand in a rolling motion around the axis of your arm to the remembered target orientation.

**Note:** In this phase the virtual tool gives no indication about the orientation of the hand around the roll axis. You must rotate the hand to the **remembered** orientation while pointing to the target in azimuth and elevation.

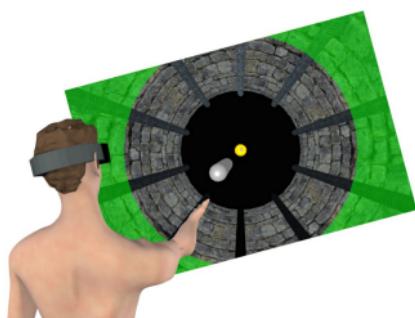


With your left hand, press the **Select Button** on the **VR Headset Remote** to validate your response once you think your hand is aligned with the memorised target.




## *Step 5: Check Response*

In some of the trials a line of balls will be projected from your hand towards the targets, allowing you to verify whether or not you correctly aligned your hand. On other trials you will see only a single spherical projectile without any indication as to whether your hand was correctly aligned to the memorised target orientation.


A green sphere appears briefly if the task sequence is completed correctly, regardless of the accuracy of the hand's alignment to the target.

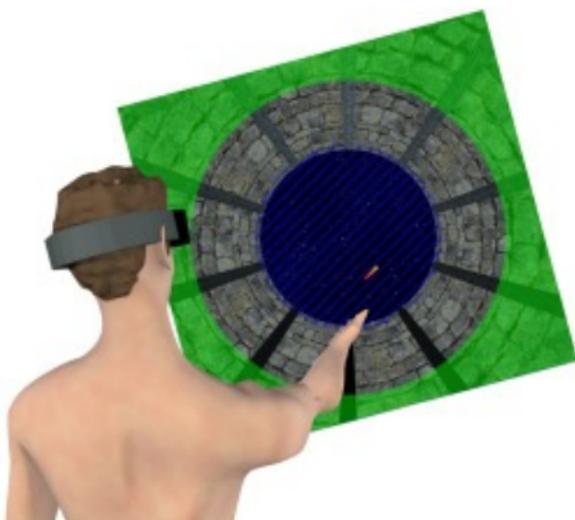

The Visual-Manual task starts with a welcome screen that prompts the subject to press a button in order to start a block of trials:

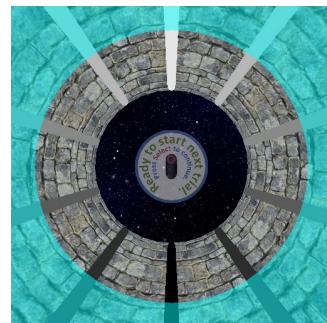


The subject is prompted to look straight ahead and bring the head to an upright position with respect to the shoulders.

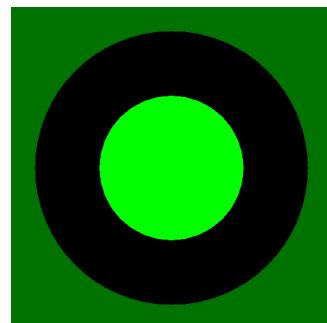


The subject should straighten the head first, then press the Select button on the VR Remote to continue.

A target orientation is then presented to the subject in the form of a grill of parallel lines situated at the end of the tunnel:

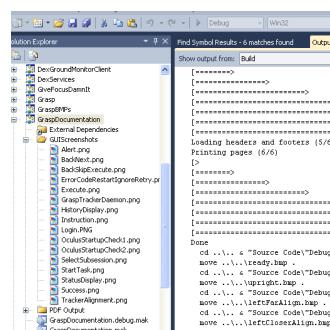


The target disappears and the subject is prompted to tilt the head to a new orientation through the same system of halo colors and direction arrows as described previously:

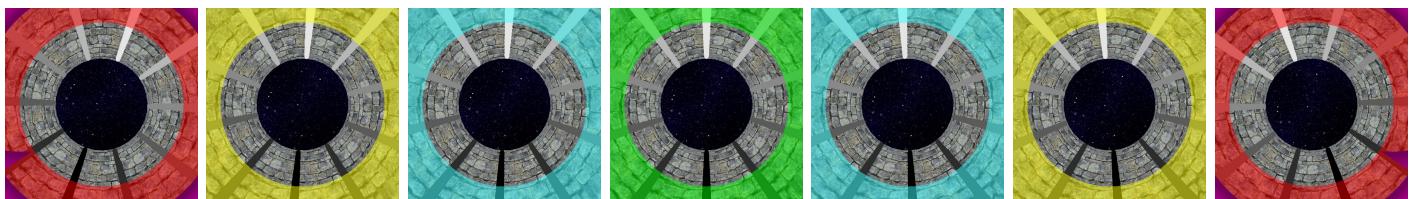

Tilted to the Left ----- Correct Orientation ----- Tilted to the Right

The subject then reproduces the remembered target orientation with the hand. The subject sees a) an orange spherical target at the end of the tunnel, b) a virtual representation of a tool in the visual scene (a cylinder) that moves with the hand in 3D visual space and c) a 'laser pointer' that projects from the tool along the axis of the fingers. The subject must raise the arm in front of the body, aim the tool by moving the arm in azimuth and elevation until the laser pointer falls on the spherical target and then rotate the palm of the hand to align with the remembered orientation of the visual target.

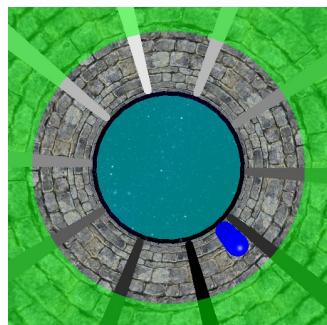


The subject presses the Select Button on the VR Remote to validate the trial and move on to the next. If the block of trials is completed, a message is displayed to let the subject know:

The screenshot shows the Visual Studio IDE interface. The Solution Explorer window on the left displays a solution named 'DevGroundMonitorClient' containing projects like 'DevServices', 'GrafFocusClient', 'GrafServer', 'GraphPPI', 'GroundMonitor', and 'GroundMonitorClient'. The Properties window in the center shows settings for the 'GroundMonitorClient' project, including build configurations like 'Debug' and 'Win32'. The Output window at the bottom right shows build logs, indicating a successful build with 0 errors and 0 warnings.

The subject may press the Select button on the VR Remote to exit the VR scene and pause the experiment. Or he or she can press and hold the Return button on the VR Remote to continue to the next block of trials.

If all blocks in the subsession are completed, the subject is prompted to doff the Virtual

## Reality Headset.

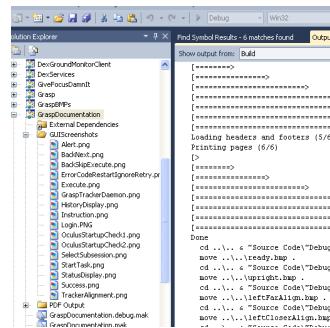

Note that other circular prompt with text instructions may appear at different times in response to different conditions or actions by the subject. For instance, a message will be generated if the subject does not maintain the head in the specified tilted position. All possible messages are shown in a separate section of this document.


## Ready to Start



Press the '**Start**' button below to continue.

## Ready to Start


- Adopt the *seated* posture in the ***GRIP Chair***.
- Attach ***Chest Marker Support*** to the ***Chair Belt Junction*** and verify that it is well centered left-to-right on your chest.
- Place the ***Virtual Reality Headset*** on your forehead in preparation for performing science tasks.
- Pass your right forearm through the strap on the ***Hand Marker Support*** cable and tighten near the elbow.
- Attach the ***Hand Marker Support*** to your ***right*** hand with your palm against the flat side of the paddle.
- Take the ***VR Headset Remote*** in your ***left*** hand.

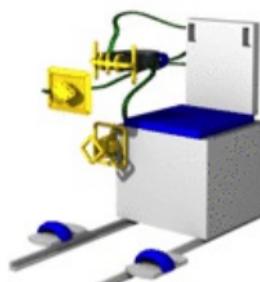

Pull the ***Virtual Reality Headset*** over your eyes then press the ***Return Button*** on the ***VR Headset Remote*** to begin.


## Ready to Start

- Adopt the *quasi-freefloating* posture in the **Floating Backrest**.
- Attach **Chest Marker Support** to your torso and verify that it is well centered left-to-right on your chest.
- Place the **Virtual Reality Headset** on your forehead in preparation for performing science tasks.
- Attach the **Hand Marker Support** to your *right* hand and attach strap on cable to forearm near elbow.
- Take the **VR Headset Remote** in your *left* hand.

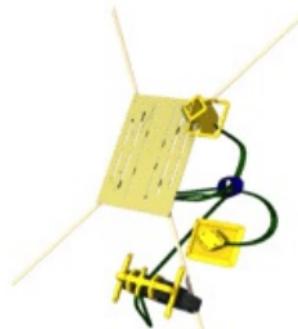

Pull the **Virtual Reality Headset** over your eyes then press the **Return Button** on the **VR Headset Remote** to begin.

## Application Running

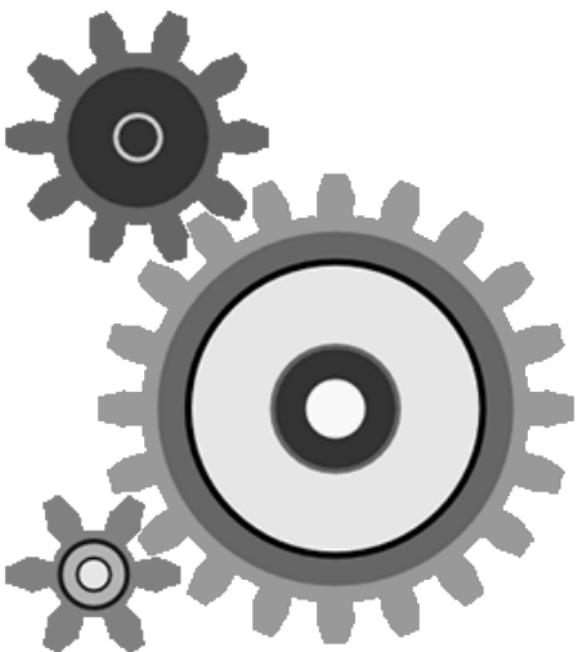

## Action Completed Normally



Press '**Next**' to continue.

## Action Completed with Errors



Press:


- '**Retry**' to repeat this action.
- '**Ignore**' to ignore and move on.
- '**Restart**' to restart this task or select another.

**Recommendation:** Attempt '**Retry**' no more than once or twice before consulting with COL-CC. '**Ignore**' and '**Restart**' should not be used without COL-CC authorization.

## Tracker Alignment

Prepare to perform the alignment procedure for the 3D tracker:

- Check that the ***GRIP Chair*** is installed with the seat back deployed and the ***Chest Marker Support, Hand Marker Support*** and ***Head Marker Support*** are attached with velcro to the ***GRIP Chair*** as shown.

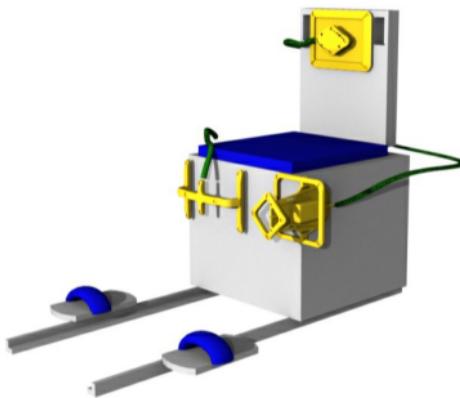

**Note:** Alignment of the ***Chest Marker Support*** on the backrest is critical. Please ensure that the top edge of the ***Chest Marker Support*** is aligned with the top edge of the backrest with the cable exiting to the right side of the chair (subject's right side when seated facing the ***Tracking Cameras***).


- Check that the field of view of the ***Tracking Cameras*** is free of obstructions.

Press '**Start**' to perform alignment.

## Tracker Alignment

Prepare to perform the alignment procedure for the 3D tracker:

- Check that the quasi-freefloating ***Floating Backrest*** is installed with the ***Chest Marker Support, Hand Marker Support*** and ***Head Marker Support*** attached with velcro as shown to the restraint plate.

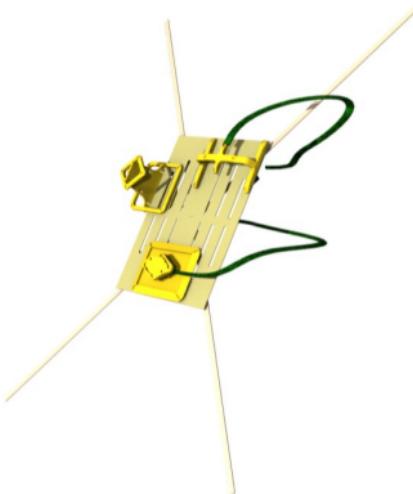

**Note:** Alignment of the ***Chest Marker Support*** on the restraint plate is critical. Check that the top edge of the ***Chest Marker Support*** is aligned with the edge of the plate nearest the deck of the ISS module and that the cable exits the support to your right when facing the plate.

- Check that the field of view of the ***Tracking Cameras*** is free of obstructions.

Press '**Start**' to perform alignment.

## Alignment Process Interrupted



The tracker alignment process was terminated prior to completion.

It **may** be acceptable to continue if a successful alignment has recently been performed **and** if the tracker cameras have not moved in the meantime.

Press:

- '**Retry**' to repeat the alignment procedure.

Or consults COL-CC and press:

- '**Ignore**' to ignore and move on.
- '**Restart**' to restart this task or select another.

**Recommendation:** Attempt '**Retry**' no more than once or twice before consulting with COL-CC. '**Ignore**' and '**Restart**' should not be used without COL-CC authorization.

# Tracker Alignment Failure



Verification of the tracker alignment failed. This may be due to occlusions in the field of view of each tracking camera, or may be due to a misconfiguration of the tracker hardware.

To continue:

- Check for possible occlusions and then press '**Retry**' to repeat this action.


If not successfull, consult with COL-CC and then press:

- '**Retry**' to repeat the alignment procedure.
- '**Ignore**' to ignore and move on.
- '**Restart**' to restart this task or select another.

**Recommendation:** Attempt '**Retry**' no more than once or twice before consulting with COL-CC. '**Ignore**' and '**Restart**' should not be used without COL-CC autorization.

## HMD Assembly

- Attach the ***Head Marker Support*** to the front of the ***Virtual Reality Headset*** as shown:

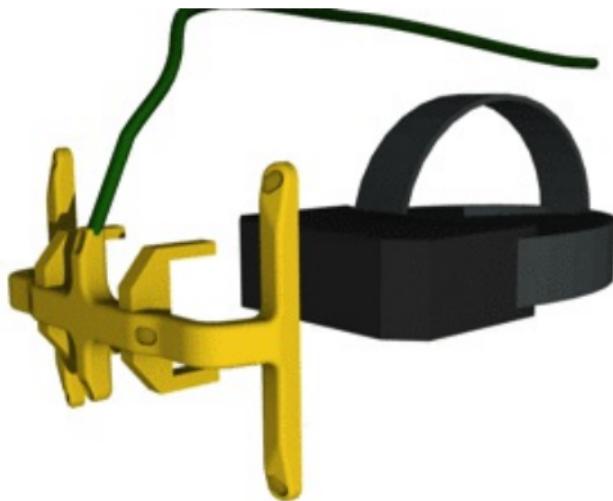

- Secure with velcro straps on top and bottom.
- Attach the small junction box on the ***Head Marker Support*** cable to the top strap of the ***Virtual Reality Headset***.


**Note:** Alignment of the Marker Structure to the ***Virtual Reality Headset*** is critical. Please ensure that the structure is parallel to the front face of the ***Virtual Reality Headset***.

Press '**Next**' to continue.

## Objects Pre-Positioning

Place the objects as shown to prepare for experiment execution:

- Replace the Cushion on the ***GRIP Chair***.
- Temp stow the ***Virtual Reality Headset***, the ***Hand Marker Support*** and the ***Chest Marker Support*** to the right side of the ***GRIP Chair*** (subject's right side when seated facing the ***Tracking Cameras***).

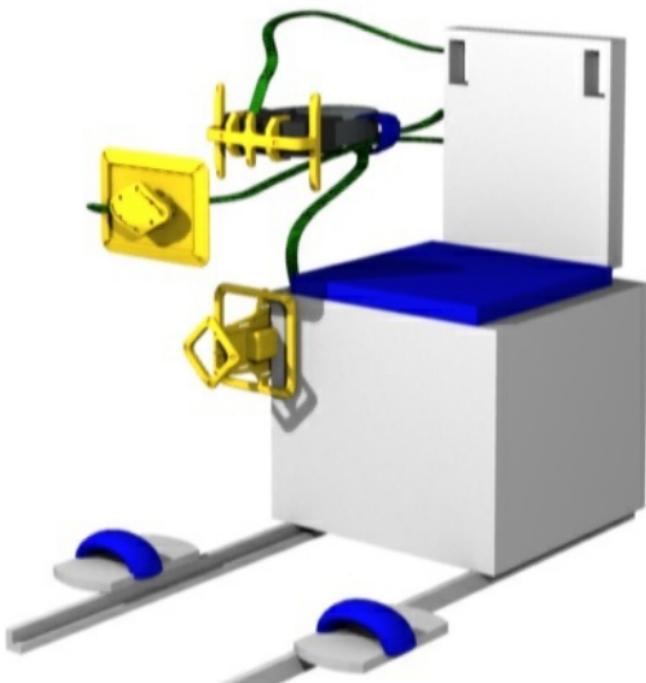


- Ensure that the ***Marker Array Cables*** and ***VR Headset Cable*** pass behind the ***GRIP Chair***.

Press 'Next' to continue.

## Object Pre-Positioning

Place the objects as shown to prepare for experiment execution:

- Temp stow the ***Virtual Reality Headset***, the ***Hand Marker Support*** and the ***Chest Marker Support*** to the right side of the ***Floating Backrest*** as shown:

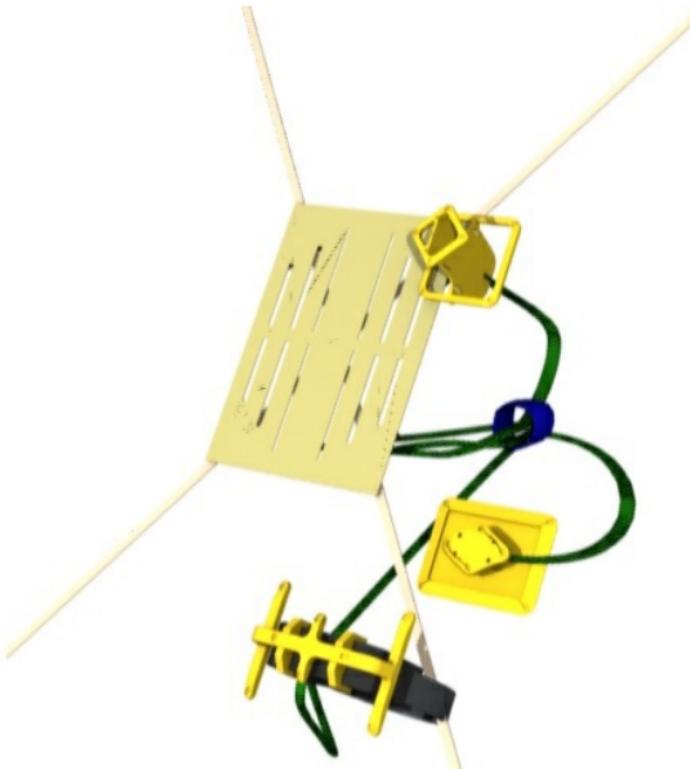

- Ensure that the ***Marker Array Cables*** and ***VR Headset Cable*** pass behind the ***Floating Backrest***.

Press 'Next' to continue.

## File Transfer Error



An error occurred when attempting to transfer files to the GRIP payload.

- Check that GRIP Chair is powered on.
- Check that Perspectives Laptop is connected to GRIP Chair.

Press:

- '**Retry**' to repeat this action.
- '**Ignore**' to ignore and move on.
- '**Restart**' to restart this task or select another.

**Recommendation:** Attempt '**Retry**' no more than once or twice before consulting with COL-CC. '**Ignore**' and '**Restart**' should not be used without Col-CC authorization.

## Incomplete Action



Execution of experiment trials was interrupted prior to completion of all trials.

Options to continue:

- Press '**Retry**' to repeat this block of trials from the beginning.
- Press '**Ignore**' to ignore and move on (per COL-CC only).
- To complete only the missing trials (per COL-CC only):
  - Note the number of the interrupted trial block (highlighted in the task list).
  - Select *Resume Interrupted* to complete the aborted trials.
  - Return to the trial block following the interrupted block in the task list and continue normally to complete the subsession.

## Max Retries Exceeded



The maximum allowable number of trials to be repeated has been exceeded. This can occur if:

- too many procedural errors (e.g. timeouts) have occurred.
- you have voluntarily aborted too many trials.
- you have voluntarily rejected too many responses.

Please consult with COL-CC and then press:

- '**Retry**' to repeat the block of trials.
- '**Ignore**' to ignore and move on.
- '**Restart**' to restart this task or select another.

**Note:** '**Ignore**' and '**Restart**' should not be used without COL-CC authorization.

## Set No GRIP Tracking

Executing this task will set a flag to tell GRASP to use only the inertial tracking inside the Oculus.

- The Visual-Visual should work normally.
- The Visual-Manual and the Manual-Manual paradigms will work in a mode where tracking of the hand is simulated with the mouse and arrow keys.

This trackerless version allows for debugging of the user interface, procedures and operations without the need for connection to the GRIP hardware.

Press '**Start**' to continue.

## Set Normal GRIP Tracking

Executing this task will reset the system to use the normal 3D tracking system for ***Virtual Reality Headset***, hand and chest. The Perspectives computer (this computer) must be connected to the GRIP hardware and the ***Virtual Reality Headset***, Hand and Chest marker arrays must be available and connected as well.

Press '**Start**' to continue.

## File Transfer

Executing this task will package and send the data files from today's sessions to GRIP for later transfer to ground.

The Perspectives computer (this computer) must be connected to the GRIP hardware and both must be running.

Press the '**Start**' button below to continue.

**NB:** Pressing 'Enter' will not work here.  
You must click on the 'Start' button below.

# Task Completed



Press '**N**ext' to proceed.

## Subsession Completed



Press '**Next**' to exit.



BlockCompleted



DemoWorking



HandRotateTimeout



HandShouldNot



HandTooSoon



HeadAlignmentTimeout



HeadMisalignment



InvalidateTrial



LowerArm



LowerArmTimeout



RaiseArm



RaiseArmTimeout



ReadyToStart



RejectTrial



ResponseTimeout



StraightenHead



TimeLimit



VRCompleted



BlockCompleted



DemoWorking



HandRotateTimeout



HandShouldNot



HandTooSoon



HeadAlignmentTimeout



HeadMisalignment



InvalidateTrial



LowerArm



LowerArmTimeout



RaiseArm



RaiseArmTimeout



ReadyToStart



RejectTrial



ResponseTimeout



StraightenHead



TimeLimit



VRCompleted